

Manufacturer: Intini Marketing Inc.
Model: Medi-Lock
ASTM Type: X

Description:

This package consists of a one-piece plastic box with a lid permanently attached by a continuous plastic hinge on one side of the lid (Figures 1 and 2).



Figure 1

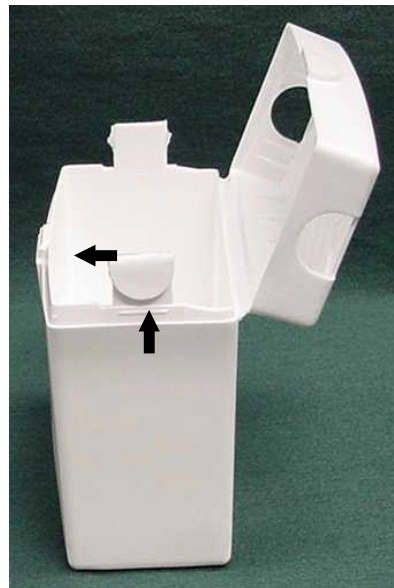


Figure 2

The box is opened by squeezing the center semicircular section on the front of the box below the lid, lifting the lid slightly until one hears a click, and then pressing the two semicircular tabs on the sides of the lid and lifting up the lid. The center semicircular section on the front of the box has a lug on the rim (see horizontal arrow in Figure 2) that fits into a semicircular indentation inside the lid when closed. The two tabs on the sides of the box each have a semicircle thicker than the tab itself (see vertical arrow in Figure 2). The two semicircles on the tabs lock into the semicircular holes on the sides of the lid when the box is closed (Figure 2).



Figure 3

The manufacturer's instructions on the lid are (Figure 3):

1. Holding base with one hand, squeeze semicircular section with your thumb.
2. With your other hand, push up cover until you hear two clicks.
3. Grasp narrow sides of the top. Using the index finger and thumb, press the two semicircular locking tabs on the sides of the cover. Lift cover completely.

(To the left of the written instructions is a series of 3 pictures showing the actions associated with the instructions.)